



### How to retube an earmould

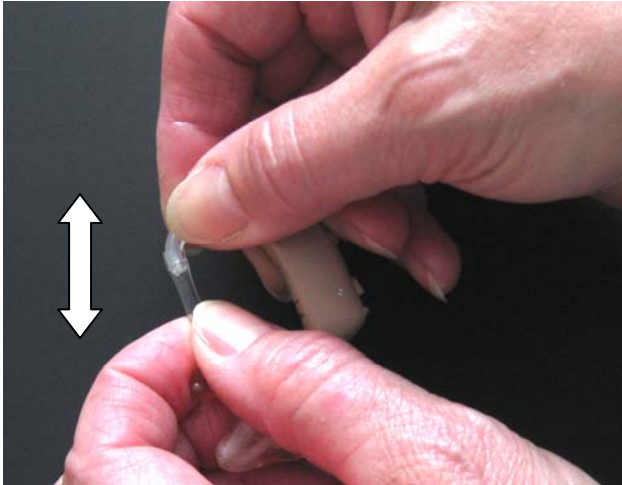

1. Pull the tube apart from the plastic elbow of the hearing aid.

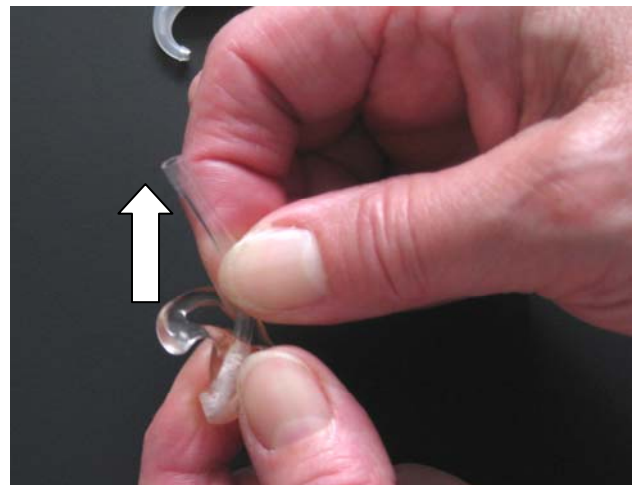

2. Pull the tube out of the earmould




3. Keep the tube to use as a measure guide for the new tube.



4. Taper the end of the new piece of tubing



5. Thread the tapered end of the tube through the hole in the earmould

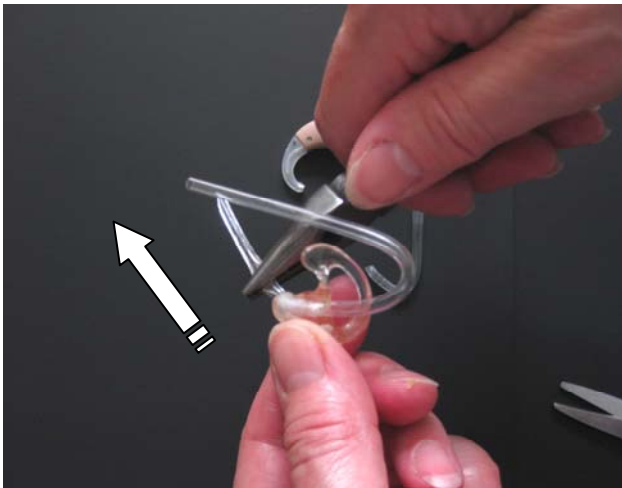

6. Carefully pull the tube through with pliers...



7. ...until the tube is in this position




8. Cut off the tapered end of the tube with a sharp blade. Make sure there are no rough edges left.



9. Use the old piece of tube as a guide as to the length that the new tube needs to be cut



10. Cut the new tube to the same length as the old tube

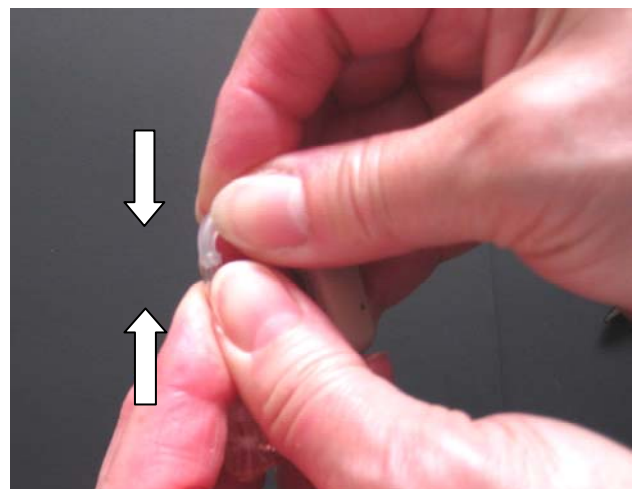

11. Push the tube over the plastic elbow of the hearing aid